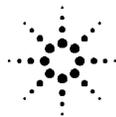

**Agilent 42030A Four-Terminal Pair Standard
Resistor Set**

Operation and Service Manual



Agilent Technologies

Innovating the HP Way

**Agilent Part No. 42030-90001
Printed in JAPAN January 2000**

Second Edition

CERTIFICATION	<p>Agilent Technologies certifies that this product met its published specifications at the time of shipment from the factory. Agilent Technologies further certifies that its calibration measurements are traceable to the United States National Institute of Standards and Technology (NIST), to the extent allowed by the Institute's calibration facility, or to the calibration facilities of other International Standards Organization members.</p>
WARRANTY	<p>This Agilent Technologies instrument product is warranted against defects in material and workmanship for a period of one year from the date of shipment. During the warranty period, Agilent Technologies will, at its option, either repair or replace products which prove to be defective.</p> <p>For warranty service or repair, this product must be returned to a service facility designated by Agilent Technologies. Buyer shall prepay shipping charges to Agilent Technologies and Agilent Technologies shall pay shipping charges to return the product to Buyer. However, Buyer shall pay all shipping charges, duties, and taxes for products returned to Agilent Technologies from another country.</p> <p>Agilent Technologies warrants that its software and firmware designated by Agilent Technologies for use with an instrument will execute its programming instruction when properly installed on that instrument. Agilent Technologies does not warrant that the operation of the instrument, or software, or firmware will be uninterrupted or error free.</p>
LIMITATION OF WARRANTY	<p>The foregoing warranty shall not apply to defects resulting from improper or inadequate maintenance by Buyer, Buyer-supplied software or interfacing, unauthorized modification or misuse, operation outside of the environmental specifications for the product, or improper site preparation or maintenance.</p> <p>NO OTHER WARRANTY IS EXPRESSED OR IMPLIED. AGILENT TECHNOLOGIES SPECIFICALLY DISCLAIMS THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE.</p>
EXCLUSIVE REMEDIES	<p>THE REMEDIES PROVIDED HEREIN ARE BUYER'S SOLE AND EXCLUSIVE REMEDIES. AGILENT TECHNOLOGIES SHALL NOT BE LIABLE FOR ANY DIRECT, INDIRECT, SPECIAL, INCIDENTAL, OR CONSEQUENTIAL DAMAGES, WHETHER BASED ON CONTRACT, TORT, OR ANY OTHER LEGAL THEORY.</p>
ASSISTANCE	<p><i>Product maintenance agreements and other customer assistance agreements are available for Agilent Technologies products.</i></p> <p><i>For any assistance, contact your nearest Agilent Technologies Sales and Service Office. Addresses are provided at the back of this manual.</i></p>

Notice The information contained in this document is subject to change without notice.

This document contains proprietary information which is protected by copyright. All rights are reserved. No part of this document may be photocopied, reproduced, or translated to another language without the prior written consent of the Agilent Technologies.

Agilent Technologies Japan, Ltd.
Component Test PGU-Kobe
1-3-2, Murotani, Nishi-ku, Kobe-shi,
Hyogo, 651-2241 Japan

Manual Printing History

The manual printing date and part number indicate its current edition. The printing date changes when a new edition is printed. (Minor corrections and updates which are incorporated at reprint do not cause the date to change.) The manual part number changes when extensive technical changes are incorporated.

October 1992 First Edition (part number: 42030-90001)

January 2000 Second Edition (part number: 42030-90001)

How To Use This Manual. This is the Operation and Service Manual for the 42030A Four-Terminal Pair Standard Resistor Set. This contains operation and service information, and consists of the following two chapters.

**Chapter 1 General
Information**

Chapter 1 describes the 42030A description, specification, and, other general information.

Chapter 2 Service

Chapter 2 describes the 42030A calibration and repair.

Contents

1. General Information

Introduction	1-1
Description	1-1
Initial Inspection	1-2
Repackaging for Shipment	1-2
Agilent Packing Materials	1-2
Non-Agilent Packing Materials	1-2
Specifications	1-3
DC Resistance	1-3
Operating Conditions	1-3
Dimensions	1-3
Weight	1-4
Supplemental Performance Characteristics	1-4
Temperature Coefficient and Stability	1-4
Storage Conditions	1-4

2. Calibration and Repair

Introduction	2-1
Calibration	2-1
Calibration Interval	2-4
Repair	2-4

Index

Tables

1-1. 42030A Contents	1-1
1-2. DC Resistance	1-3
1-3. Operating Conditions	1-3
1-4. Dimensions	1-3
1-5. Weight	1-4
1-6. DC Resistance Temperature Coefficient and Stability	1-4
1-7. Storage Conditions	1-4
2-1. 42030A Calibration Parameters	2-1
2-2. 42030A DC Resistance Calibration Uncertainty	2-2
2-3. 42035A 10 Ω Resistor AC Calibration Uncertainty	2-2
2-4. 42036A 100 Ω Resistor AC Calibration Uncertainty	2-3
2-5. 42037A 1 k Ω Resistor AC Calibration Uncertainty	2-3
2-6. 42038A 10 k Ω Resistor AC Calibration Uncertainty	2-3
2-7. 42039A 100 k Ω Resistor AC Calibration Uncertainty	2-4

General Information

Introduction

This chapter provides the 42030A Four-Terminal Pair Standard Resistor Set description, specifications, and related general information.

Description

The 42030A Four-Terminal Pair Standard Resistor Set consists of nine precision resistor standards. The set is constructed using high stability resistors and are virtually unaffected by changes in ambient temperature, ensuring measurement repeatability and reliability. Table 1-1 lists the 42030A contents.

Table 1-1. 42030A Contents

Model No. or Agilent Part No.	Description	Qty.
42031A	1 m Ω Standard Resistor	1
42032A	10 m Ω Standard Resistor	1
42033A	100 m Ω Standard Resistor	1
42034A	1 Ω Standard Resistor	1
42035A	10 Ω Standard Resistor	1
42036A	100 Ω Standard Resistor	1
42037A	1 k Ω Standard Resistor	1
42038A	10 k Ω Standard Resistor	1
42039A	100 k Ω Standard Resistor	1
42030-60100	Carrying Case	1
42030-90001	Operation and Service Manual	1
Not Assigned	Calibration Report	1

The 42030A is designed for use in calibrating the Agilent LCR meters which have four-terminal pair measurement terminals, and can be directly connected to the measurement terminals. The lcr meter documentation will give detailed information on using the 42030A.

Initial Inspection

Inspect the shipping container for damage. If the shipping container or cushioning material is damaged, it should be kept until the shipping contents have been checked for completeness and the instrument has been checked mechanically and electrically. The shipping contents should consist of the items shown in Table 1-1. If the shipment is incomplete, or if there is mechanical damage or other defects, notify your nearest Agilent Technologies office. If the shipping container is damaged, or the cushioning material shows signs of stress, notify the carrier as well as the Agilent Technologies office. Keep the shipping materials for the carrier's inspection. Your Agilent Technologies office will arrange for repair or replacement, without waiting for the claim settlement.

Repackaging for Shipment

Agilent Packing Materials

Containers and materials identical to those used in factory packaging are available from Agilent Technologies. If the instrument is being returned to Agilent Technologies for servicing, attach a tag indicating the type of service required, return address, model number, and full serial number. Also, mark the container FRAGILE to ensure careful handling. In any correspondence, refer to the unit by model number and full serial number.

Non-Agilent Packing Materials

The following general instructions should be used for re-packing with commercially available packing materials:

- a. If shipping to Agilent Technologies office or service center, attach a tag indicating the type of service required, return address, model number, and full serial number.
- b. Use a strong shipping container. A double-walled carton made of 350 pound test material is adequate.
- c. Use enough shock absorbing material (3- to 4-inch layer) around all sides of the case to provide a firm cushion and to prevent movement inside the container.
- d. Seal the shipping container securely.
- e. Mark the shipping container FRAGILE to ensure careful handling.
- f. In any correspondence, refer to the unit by model number and full serial number.

Specifications

This section gives the complete specifications for the 42030A Four-Terminal Pair Standard Resistor Set. When the 42030A is shipped from the factory, it meets the specifications listed in this section.

DC Resistance

Table 1-2 lists the nominal DC resistance for each resistance module of the 42030A.

Table 1-2. DC Resistance

Model Number	DC Resistance ¹
42031A	1 mΩ ± 0.2%
42032A	10 mΩ ± 0.2%
42033A	100 mΩ ± 0.2%
42034A	1 Ω ± 0.2%
42035A	10 Ω ± 0.1%
42036A	100 Ω ± 0.1%
42037A	1 kΩ ± 0.1%
42038A	10 kΩ ± 0.1%
42039A	100 kΩ ± 0.1%

¹ Specified under ambient environmental conditions of: Temperature 23 °C ± 5 °C , Relative Humidity ≤ 70%

Operating Conditions

The 42030A must be operated under the ambient environmental conditions listed in Table 1-3

Table 1-3. Operating Conditions

Temperature	23 °C ± 5 °C
Relative Humidity	≤ 70%

Dimensions

Table 1-4 shows the 42030A dimensions.

Table 1-4. Dimensions

Resistance Module	94 mm (W) × 31 mm (H) × 67 mm (D)
Carrying Case	400 mm (W) × 120 mm (H) × 290 mm (D)

Weight Table 1-5 shows the 42030A weight.

Table 1-5. Weight

Resistance Module	170 g (ea.)
42030A Total	3,700 g

Supplemental Performance Characteristics

This section gives the supplemental performance characteristics for the 42030A Four-Terminal Pair Standard Resistor Set . These supplemental performance characteristics are not guaranteed, they are the 42030A typical characteristics that may be expected.

Temperature Coefficient and Stability

Table 1-6 shows the 42030A DC resistance temperature coefficient and stability.

**Table 1-6.
DC Resistance Temperature Coefficient and Stability**

Model Number	Nominal Value	Temperature Coefficient¹	Stability²
42031A	1 m Ω	± 15 ppm/ $^{\circ}$ C	± 100 ppm/year
42032A	10 m Ω	± 15 ppm/ $^{\circ}$ C	± 100 ppm/year
42033A	100 m Ω	± 15 ppm/ $^{\circ}$ C	± 100 ppm/year
42034A	1 Ω	± 15 ppm/ $^{\circ}$ C	± 100 ppm/year
42035A	10 Ω	± 10 ppm/ $^{\circ}$ C	± 100 ppm/year
42036A	100 Ω	± 10 ppm/ $^{\circ}$ C	± 100 ppm/year
42037A	1 k Ω	± 10 ppm/ $^{\circ}$ C	± 100 ppm/year
42038A	10 k Ω	± 10 ppm/ $^{\circ}$ C	± 100 ppm/year
42039A	100 k Ω	± 10 ppm/ $^{\circ}$ C	± 100 ppm/year

¹ Under temperature condition: 23 $^{\circ}$ C \pm 5 $^{\circ}$ C

² Under environmental conditions: Temperature 23 $^{\circ}$ C \pm 5 $^{\circ}$ C , Relative Humidity \leq 70%

Storage Conditions

The 42030A must be stored or shipped under the ambient environmental conditions listed in Table 1-7.

Table 1-7. Storage Conditions

Temperature	- 40 $^{\circ}$ C to + 70 $^{\circ}$ C
Relative Humidity	\leq 95% at 40 $^{\circ}$ C

Calibration and Repair

Introduction

This chapter provides the calibration and repair information for the 42030A Four-Terminal Pair Standard Resistor Set .

Calibration

The 42030A is calibrated in the parameters given in Table 2-1 with the calibration uncertainties given in Table 2-2 through Table 2-7 and equipped with a calibration report, when shipped from the factory.

Table 2-1. 42030A Calibration Parameters

Model Number	Nominal Value	DC Resistance	Series Resistance and Reactance	Parallel Resistance and Susceptance
42031A	1 m Ω	✓		
42032A	10 m Ω	✓		
42033A	100 m Ω	✓		
42034A	1 Ω	✓		
42035A	10 Ω	✓	✓	
42036A	100 Ω	✓	✓	
42037A	1 k Ω	✓		✓
42038A	10 k Ω	✓		✓
42039A	100 k Ω	✓		✓

Table 2-2.
42030A DC Resistance Calibration Uncertainty

Model Number	Nominal Value	Calibration Uncertainty
42031A	1 m Ω	$\pm 0.1\%$
42032A	10 m Ω	$\pm 0.1\%$
42033A	100 m Ω	$\pm 0.1\%$
42034A	1 Ω	$\pm 0.1\%$
42035A	10 Ω	$\pm 0.03\%$
42036A	100 Ω	$\pm 0.02\%$
42037A	1 k Ω	$\pm 0.02\%$
42038A	10 k Ω	$\pm 0.02\%$
42039A	100 k Ω	$\pm 0.02\%$

Table 2-3.
42035A 10 Ω Resistor AC Calibration Uncertainty

Frequency	Calibration Uncertainty	
	Series Resistance	Reactance
1 MHz	$\pm 0.12\%$	$\pm 0.010 \Omega$
2 MHz	$\pm 0.15\%$	$\pm 0.010 \Omega$
3 MHz	$\pm 0.20\%$	$\pm 0.015 \Omega$
4 MHz	$\pm 0.25\%$	$\pm 0.025 \Omega$
5 MHz	$\pm 0.4\%$	$\pm 0.03 \Omega$
10 MHz	$\pm 1.0\%$	$\pm 0.08 \Omega$
13 MHz	$\pm 1.5\%$	$\pm 0.10 \Omega$

Table 2-4.
42036A 100 Ω Resistor AC Calibration Uncertainty

Frequency	Calibration Uncertainty	
	Series Resistance	Reactance
1 MHz	$\pm 0.12\%$	$\pm 0.10 \Omega$
2 MHz	$\pm 0.12\%$	$\pm 0.10 \Omega$
3 MHz	$\pm 0.15\%$	$\pm 0.15 \Omega$
4 MHz	$\pm 0.15\%$	$\pm 0.15 \Omega$
5 MHz	$\pm 0.20\%$	$\pm 0.20 \Omega$
10 MHz	$\pm 0.5\%$	$\pm 0.5 \Omega$
13 MHz	$\pm 0.7\%$	$\pm 0.7 \Omega$

Table 2-5.
42037A 1 k Ω Resistor AC Calibration Uncertainty

Frequency	Calibration Uncertainty	
	Parallel Resistance	Susceptance
100 kHz	$\pm 0.12\%$	$\pm 0.8 \mu\text{S}$
1 MHz	$\pm 0.08\%$	$\pm 1.0 \mu\text{S}$
2 MHz	$\pm 0.10\%$	$\pm 1.0 \mu\text{S}$
3 MHz	$\pm 0.10\%$	$\pm 1.5 \mu\text{S}$
4 MHz	$\pm 0.15\%$	$\pm 2.0 \mu\text{S}$
5 MHz	$\pm 0.15\%$	$\pm 3.0 \mu\text{S}$
10 MHz	$\pm 0.4\%$	$\pm 5 \mu\text{S}$
13 MHz	$\pm 0.6\%$	$\pm 8 \mu\text{S}$

Table 2-6.
42038A 10 k Ω Resistor AC Calibration Uncertainty

Frequency	Calibration Uncertainty	
	Parallel Resistance	Susceptance
100 kHz	$\pm 0.08\%$	$\pm 0.10 \mu\text{S}$
1 MHz	$\pm 0.08\%$	$\pm 0.20 \mu\text{S}$

Table 2-7.
42039A 100 k Ω Resistor AC Calibration Uncertainty

Frequency	Calibration Uncertainty	
	Parallel Resistance	Susceptance
100 kHz	$\pm 0.08\%$	$\pm 0.02 \mu S$
1 MHz	$\pm 0.12\%$	$\pm 0.15 \mu S$

The 42030A must be calibrated at a facility which satisfies the calibration uncertainties given in Table 2-2 through Table 2-7. Agilent Technologies's calibration laboratories satisfy all of these calibration uncertainties. For complete information on calibration, contact your nearest Agilent Technologies sales and service office.

Calibration Interval The 42030A's recommended calibration interval is one year.

Repair

The 42030A can be repaired only at Agilent Technologies facilities. For complete information on repair, contact your nearest Agilent Technologies Sales and Service Office.